



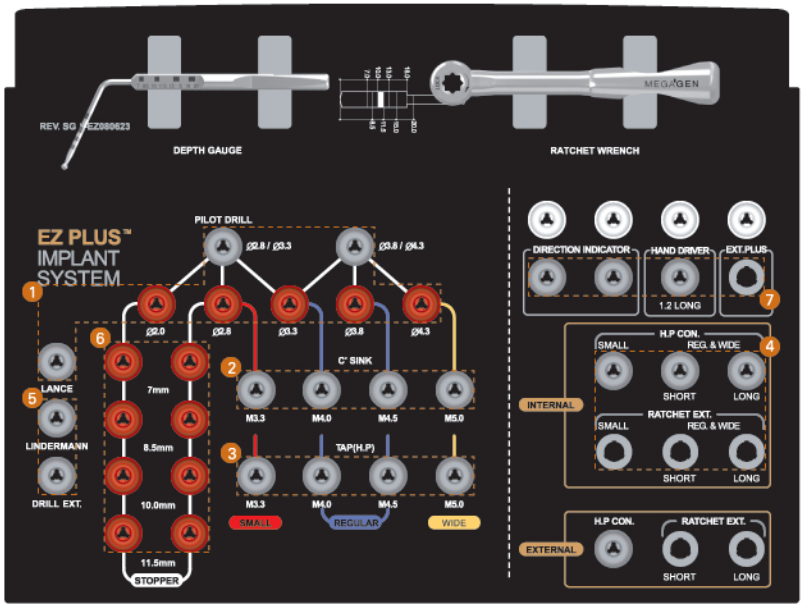

[FDA clearance # K070562]

MegaGen USA Inc. [NJ] Tel: 1-201-678-2457 / Toll Free: 866-277-5662 / E-mail: megagenusa@gmail.com MegaGen Implant Inc. [CA] Tel: 1-213-480-5770 / E-mail: megageninc@gmail.com
MegaGen Implant Hong Kong Limited Tel: 852-2866-0199 / E-mail: mynameisjwon@gmail.com
MegaGen Implant Singapore PTE. LTD Tel: 65-6519-3359 / E-mail: hskim1978@gmail.com
MegaGen Kyrgyzstan Tel: 996-312-646017 / E-mail: wametro@naver.com MegaGen Tashken Tel: 998-71-256-2720 / E-mail: yugcho@hanmail.net
MegaGen Implant Co.,Ltd [Head Office] 377-2, Kyochon-Ri, Jain-Myun, Gyeongsbuk, Gyeongsbuk, Korea, 712-852 / Tel: 82-53-857-5770 / E-mail: sales@megagen.co.kr

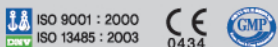

EZ Plus™ Internal

Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen!

Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen!



MEGA^{GEN} IMPLANT

Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen!

Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen! Dental Implant Special Company,MegaGen!

MEGAGEN



EZ Plus™ Internal

Get the Maximum Benefits

Submerged Type Implant, Trip Internal Connection

Design The EZ Plus fixture's crestal tapered design provides maximum initial stability in any bone quality. In dense bone situations, the self tapping thread design and four cutting edges of our fixture provide easy insertion. The thread design that runs to the top of the fixture has been engineered to maximize the connection between the bone and the implant, while minimizing crestal bone loss. Platform switching enhances crestal bone attachment.

Surface Our Super RBM Surface Treatment produces a favorable surface topography that has been proven through extensive systematic testing. MegaGen's Super RBM Surface shortens loading time and provides favorable results in any bone quality.


Feature Internal Trip Connection : MegaGen's Trip anti-rotation structure will not deform until 220N of insertion torque, allowing for extremely stable initial fixation.

Surgery On a case by case basis, proper initial fixation can be achieved using the EZ Plus Counter Sink Drill. The mount-free type handpiece connector with a built-in ceramic ball provides stability so that fixture adhesion can be maintained even after many repeated uses.

Prosthetics The EZ Plus system incorporates an 11° Morse taper, producing a conical seal which forms a cold weld between the abutment and the implant. Abutment options to provide esthetic and prosthetic excellence include UCLA, Angled, Solid, Post, Milling type, Ceramic, and O-Ring abutments as well as various other abutment options.

Packaging

Pre-Mount Type : EZ Plus Internal Fixture Connected with Multi-Mount and Cover Screw
No-Mount Type : EZ Plus Internal Fixture and Cover Screw

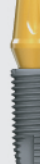	Small	Regular				Wide	
							
Type	MF	MF	MB	MF	MB	MF	MB
Fixture Diameter	Ø3.3	Ø4.0		Ø4.5		Ø5.0	
Body Length		7.0 / 8.5 / 10 / 11.5 / 13 / 15mm					

Strong Internal Trip™ (Triple Groove)



Convenience of Mount Free

EZ Plus MF is a Mount Free type implant. Free from inconvenience to remove the mount after placing the implant. Enjoy the convenience of Mount Free & Economic Price.

